11 Drugs Rapid Test Panel (Urine) Package Insert For Self-testing

A rapid test for the qualitative detection of multiple drugs and drug metabolites in human urine.
For self-testing, not for in vitro diagnostic use.

【INTENDED USE】

The 11-Drug Rapid Test Panel is a rapid chromatographic immunoassay for the qualitative detection of multiple drugs and drug metabolites in urine at the following cut-off concentrations:

Test	Calibrator	Cut-off (ng/mL)
Amphetamine (AMP)	d-Amphetamine	1,000
Barbiturates (BAR)	Secobarbital	300
Benzodiazepines (BZO)	Oxazepam	300
Cocaine (COC)	Benzoylecgonine	300
Marijuana (THC)	11-nor-Δ⁹-THC-9 COOH	50
Methadone (MTD)	Methadone	300
Methamphetamine (MET)	d-Methamphetamine	1,000
Methylenedioxy-methamphetamine (MDMA)	d,l-Methylenedioxymethamphetamine	500
Morphine/Opiate (MOP/OPI)	Morphine	300
Fentanyl (FYL)	Fentanyl	10
Mephedrone (MEP)	Mephedrone	100

The test provides only preliminary test results. A more specific alternative chemical method must be used in order to obtain a confirmed analytical result. GC/MS or LC/MS is the preferred confirmatory method.

【SUMMARY】

The 11-Drug Rapid Test Panel is a rapid urine screening test that can be performed without the use of an instrument. The test utilizes monoclonal antibodies to selectively detect elevated levels of specific drugs in urine.

【PRECAUTIONS】

- Read the entire package insert prior to performing test.
- For self-testing, not for *in vitro* diagnostic use.
- For external use only.
- Do not use the test after expiration date printed on the package.
- Do not use the test if its foil pouch is torn or damaged.
- For single use. Discard after first use.
- The test device should remain in the sealed pouch until use.
- Contaminated or tainted urine sample may give false results.
- Keep out of reach of children.

【MATERIALS PROVIDED】

- Test Panel (sealed in foil pouch with a desiccant)
- Package Insert

【MATERIALS REQUIRED BUT NOT PROVIDED】

- Watch, clock or timer
- Urine collection cup

【STORAGE AND STABILITY】

Store as packaged in the sealed pouch at 2-30 °C up to the date of expiration.
Keep away from direct sunlight, moisture and heat. DO NOT FREEZE.

【DIRECTIONS FOR USE】

Allow the test to reach room temperature (15-30 °C) before testing.
Wash your hands thoroughly.

Step 1: Collect Urine Sample

- Collect urine in a clean and dry container.


Step 2: Testing

- Remove the test panel from the sealed pouch and use it immediately.
- Remove the cap of the panel.
- Immerse the test panel vertically in the urine specimen for at least **10 to 15 seconds**.
- **Note:** Immerse the dipstick panel to at least the level of the wavy lines, but not above the arrow on the test panel.
- Replace the cap and place the test panel on a non-absorbent flat surface.
- Start a timer.


Step 3: Read Results

- **Read the test results at 5 minutes,** do not read results after 10 minutes.

- Remove the cap.
- Immerse 10-15 seconds

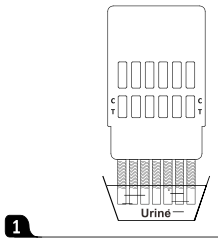1

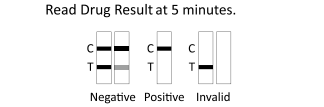2

【INTERPRETATION OF RESULTS】

	**NEGATIVE: Two colored lines appear.** One line appears next to the letter 'C' and the other line appears next to the letter 'T'.
	**POSITIVE: Only one colored line appears next to the letter 'C'. No colored line appears next to the letter 'T'.**
	**INVALID: Control line fails to appear.**

IMPORTANT: The result you obtained is called preliminary for a reason. The sample must be tested by laboratory in order to confirm the result.

【Q&A】

1. How do I know if the test worked well?
When the control line(C) appears, it means that the test unit is working well.
2. How soon can I read my results?
You can read your results at 5 minutes as long as a colored line has appeared next to the Control region(C), do not read results after 10 minutes.
3. When is the best time to run the test?
You can run the test at any time of the day.
4. How to read the test if the color and the intensity of the lines are different?
The color and intensity of the lines have no importance for result interpretation. The test should be considered as negative whatever the color intensity of the test line (T) is.
5. What is a false positive test?
A false positive test result means the drug is not present but shows detected by the device. The most common causes of a false positive test are cross reactants.
6. What is a false negative test?
A false negative test means the drug is present but is not detected by the device. If the sample is diluted, or the sample is contaminated that may cause a false negative result.

【LIMITATIONS】

1. The test is for use with human urine only.
2. The test only provides a preliminary result.
3. Substances, such as bleach and/or alum, in urine specimens may produce incorrect results. If adulteration is suspected, the test should be repeated with another urine specimen.
4. Test does not distinguish between drugs of abuse and certain medications.
5. A positive test result may be obtained from certain foods or food supplements.

Index of Symbols

	Consult instructions for use or consult electronic instructions for use		Contains sufficient for <n> tests	30°C 2°C	Temperature limit
	Caution	LOT	Batch code	REF	Catalogue number
EU REP	European Authorized Representative		Use-by date		Do not re-use
	Do not use if package is damaged and consult instructions for use		Manufacturer		

Hangzhou AllTest Biotech Co.,Ltd.
#550,Yinhai Street
Hangzhou Economic & Technological Development Area
Hangzhou, 310018 P.R. China
Web: www.alltests.com.cn Email: info@alltests.com.cn

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MedNet EC-REP GmbH